


Separation and Characterization of Nanoplastics in the Environment

publication date: Sep 5, 2019 | author/source: Postnova Analytics GmbH

[Previous](#) | [Next](#)

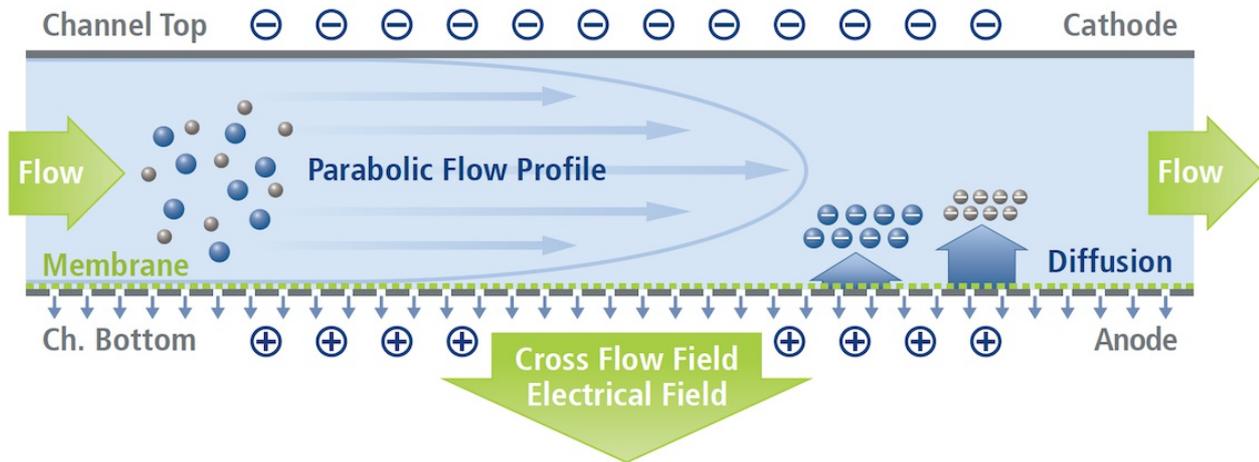

Postnova Analytics has published a new application note that demonstrates the promise of Electrical Asymmetric Field Flow Fractionation (EAF4) technology for separation and characterization of nanoplastics in the environment.

The widespread use of plastics has created one of the worst human-created environmental disasters of current times.

Typically formed by the weathering and breakdown of plastic materials, small fragments of plastic smaller than a one micrometre in size (nanoplastics) have been found to be widespread in the environment. Nanoplastics are very difficult to isolate from their environment with simple methods, such as filtration, that can be used for larger size particles. Being sub-micrometre in size, nanoplastic particles

can potentially penetrate tissues and therefore pose a significant environmental threat to marine life. However, many questions still need to be answered to establish just how dangerous nanoplastic particles really are, and a reliable method of analysis is a key requirement.



Nanoplastic particles have proven challenging to separate and characterize using traditional techniques such as dynamic light scattering or size exclusion chromatography. By comparison - Asymmetrical Flow Field-Flow Fractionation (AF4) may be the analytical solution for nanoplastics as it is a high-resolution separation technique proven for characterization of both monodisperse and polydisperse nanoparticle suspensions.

The authors of the new application describe how a mixture of two polystyrene particles was used as a test standard for a polydisperse nanoplastics system. This mixture was separated using the AF4 technique with an additional electrical field applied across the separation channel which separates the particles based on their charge polarity and magnitude. By using four different electrical field conditions, measurement of the electrophoretic mobility and thus the individual surface zeta potential of both particles in the

mixture was determined. In addition, Multi Angle Light Scattering (MALS) was used as a detector to simultaneously collect information about the size of both particles.

Compared to traditional 'batch' analysis, where only a single average size and charge is given, this new Electrical Asymmetric Field Flow Fractionation (EAF4) technique offers great promise as a powerful tool for the accurate characterisation of nanoplastics in the environment.

To request a copy of the application note please contact Postnova Analytics on +49-8191-985-6880 / +44-1885-475007 / +1-801-521-2004.

[Principle and Advantages of EAF4](#)

About Postnova Analytics

Founded in 1997, Postnova Analytics is the inventor and leading international supplier of Field-Flow Fractionation (FFF) systems for markets including biopharmaceuticals, polymers, materials, nanotechnology and environmental sciences. Leveraging its unique and patented modular FFF - Light Scattering Platform, Postnova has been able to provide high performance solutions to a wide range of applications. Today the company's growing product portfolio includes Flow FFF, Centrifugal FFF, Thermal FFF, Size Exclusion Chromatography (SEC), MALS and DLS Light Scattering. Postnova is located and headquartered in Landsberg am Lech (Germany) and has subsidiaries in the USA, UK and Finland. These offices, in conjunction with a highly qualified and trained distribution network, provide informed local support to customers worldwide.

[View Company Profile and Recent News](#)

Already a member? [Log in](#)

Interested? Require further information?

Note. Your details will be referred to the company and they will provide you with more information regarding your enquiry directly

If you have not logged into the website then please enter your details below.

About You

Prof Dr Mr Mrs Miss Ms

Send Information To

I am looking to purchase this product in:



Please upgrade to a [supported browser](#) to get a reCAPTCHA challenge.

[Why is this happening to me?](#)

Related Articles:

- [New Separation Technique Solves Complex Applications](#)
- [Determining the Physicochemical and Biophysical Attributes of Monoclonal Antibodies](#)
- [High Performance MALS Detector Enables Advanced Characterization](#)
- [Comparative Analysis of a Complex Monoclonal Antibody](#)
- [Quantifying Adeno-Associated Virus Aggregation](#)
- [Fast High Throughput Protein Purification](#)
- [Field Flow Fractionation Improves Characterization of Nanoformulated Medicines](#)
- [New Method for Characterising Liposome Formulated Drugs](#)
- [Thermal Field Flow Fractionation Copolymer Case Study](#)
- [Analysis of Nanoparticles in Commercial Sunscreens](#)
- [EAF4 Technique Helps Reveal Attributes of Monoclonal Antibodies](#)
- [Fast Measurement of Nanocrystal Length & Thickness](#)
- [Separation & Characterisation of Self Assembling Block Copolymers](#)
- [Postnova Analytics Announces 2018 Seminar on Advanced Separation Techniques](#)
- [EAF2000 Advances Complex Charged Molecule Studies](#)
- [Combination Technique Improves Polysaccharide Characterization](#)
- [Characterising Nanomaterials in Food by FFF-ICP-MS Postnova Analytics](#)
- [Unlocking the Puzzle of Neurodegenerative Diseases](#)
- [MALS Detector for Size Exclusion Chromatography and Field Flow Fractionation](#)

Newsletter Sign up

Subscribe here

Subscribe to receive our newsletters for the latest news on new laboratory products, research, Industry news and more



Weekly Update | Separation Science | Microscopy & Image Analysis | Monthly Update

Popular this Month...

Our Top 10 most popular articles this month

Today's Picks...

Looking for a Supplier?

Search by company or by product

Company Name:

Product:

S E A R C H

Please note Lab Bulletin does not sell, supply any of the products featured on this website. If you have an enquiry, please use the contact form below the article or company profile and we will send your request to the supplier so that they can contact you directly.

Lab Bulletin is published by newleaf marketing communications ltd

[Previous](#) | [Next](#)

[Back to top](#)